For Research Use Only

## TIMP1 Recombinant Matched Antibody Pair, PBS Only

www.ptgcn.com

Catalog Number: MP00245-4

**Capture Antibody** Information

Catalog Number: 83151-3-PBS Host:

Rabbit human GenBank: Isotype BC000866

**Purification Method:** Protein A purification Conjugate: Unconjugated Full name:

TIMP metallopeptidase inhibitor 1

7076 Affinity:  $K_D = 2.16 \times 10^{-11}$ 

Gene ID:

 $K_{Off} = 4.92 \times 10^{-5}$  $K_{On} = 2.28 \times 10^6$ 

Conjugate:

7076

Unconjugated

**Detection Antibody** Information

Catalog Number: 83151-7-PBS Rabbit

human Isotype: GenBank: IgG

Reactivity: TIMP metallopeptidase inhibitor 1

BC000866

**Applications** 

**Tested Applications:** Sandwich ELISA

**Purification Method:** Protein A purification

Clone ID:

230329E8

Reactivity:

Clone ID:

230329F6

15.6-1000 pg/mL (Sandwich ELISA)

Recommended Dilutions:

It is recommended that this reagent should be titrated in each testing system to obtain optimal results.

**Product Information** 

MP00245-4 targets TIMP1 in immunoassays as a matched antibody pair. Validated in Sandwich ELISA.

Capture antibody: TIMP1 Recombinant antibody, PBS Only (Capture) 83151-3-PBS (230329E8). 100 µg. Concentration 1 mgl/ml.

Detection antibody: TIMP1 Recombinant antibody, PBS Only (Detector) 83151-7-PBS (230329F6). 100 μg. Concentration 1 mgl/ml.

Alternative TIMP1 matched antibody pairs: MP00245-1, MP00245-2, MP00245-3, MP50084-1, MP50084-2

 $Unconjugated\ rabbit\ recombinant\ monoclonal\ antibody\ pair\ in\ PBS\ only\ storage\ buffer\ at\ a\ concentration\ of\ 1$ mg/mL, ready for conjugation. Created using Proteintech's proprietary in-house recombinant technology. Recombinant production enables unrivalled batch-to-batch consistency, easy scale-up, and future security of supply.

Matched antibody pairs are designed for use in a variety of assays and platforms that require matched antibody

Antibody use should be optimized for each application and assay.

Storage


Storage:


Store at -80°C.


The product is shipped with ice packs. Upon receipt, store it immediately at -80°C

PBS only


## **Selected Validation Data**

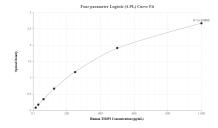

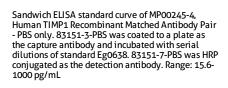

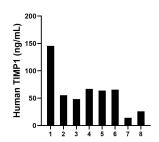

Saliva of eight individual healthy human donors was measured. The human TIMP1 concentration of detected samples was determined to be 60.8 ng/mL with a range of 14.2 - 145.8 pg/mL

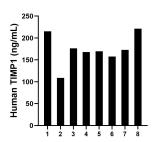

Serum of eight individual healthy human donors was measured. The human TIMP1 concentration of detected samples was determined to be 173.6 ng/mL with a range of 108.8 - 221.2 pg/mL